Operating Instructions | Betriebsanleitung | Mode d'emploi | Manuale d'uso | Instrucciones de manejo | Руководство по эксплуатации | 操作说明 | 取扱説明書

Original Operating Instructions | Original-Betriebsanleitung | Mode d'emploi original | Manuale d'uso originale | Instrucciones de manejo originales | 原始操作说明 | 取扱説明書 (オリジナル版)

YDP30 | YDP30-NET

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Contents

1 About this Document	5	6 System Settings	15
1.1 Scope.....	5	6.1 Configure "RS232" Connection (YDP30 Only).....	15
1.2 Target Groups.....	5	7 Cleaning and Maintenance	16
1.3 Symbols Used.....	5	7.1 Cleaning Device Surfaces.....	16
1.3.1 Warnings in Operation Descriptions.....	5	7.2 Maintenance Schedule.....	16
1.3.2 Other Symbols.....	5	7.3 Removing the Ink Ribbon.....	16
2 Safety Instructions	5	7.4 Removing Media from Media Mounts.....	17
2.1 Intended Use.....	5	7.5 Cleaning the Sensors and Interior Components.....	17
2.1.1 Modifications to the Device.....	6	7.5.1 Cleaning the Print Head.....	17
2.1.2 Device Repairs and Maintenance.....	6	7.5.2 Cleaning the Print Roller or Tear Bar.....	17
2.2 Personnel Qualification.....	6	7.5.3 Cleaning Sensors or the Interior.....	17
2.3 Significance of these Instructions.....	6	7.6 Calibrating Sensors and Changing Settings.....	18
2.4 Functionality of the Device.....	6	7.6.1 Check the Configuration Data	
2.5 Electrical Equipment.....	6	of the Device.....	19
2.5.1 Damage to the Device's Electrical Equipment.....	6	7.6.2 Checking Core Dump.....	19
2.5.2 Power Supply Unit and Power Supply Cable.....	6	7.7 Re-starting the System.....	19
2.6 Conduct in an Emergency.....	7	8 Malfunctions	20
2.7 Accessories, Consumables, and Spare Parts.....	7	8.1 Device Malfunctions.....	20
3 Device Description	8	8.2 Malfunctions when Printing.....	20
3.1 Device Overview.....	8	8.3 Poor Print Quality.....	21
3.2 Connections.....	8	8.4 Re-starting the System.....	21
3.3 LED Display.....	8	9 Decommissioning	21
3.4 Components for Print Function and		9.1 Decommissioning the Device.....	21
Lock for Top Cover.....	9	10 Transport	22
3.5 Components for Media Feed		10.1 Transporting the Device.....	22
and Sensors.....	9	11 Storage and Shipping	22
4 Installation	10	11.1 Storing.....	22
4.1 Scope of Delivery.....	10	11.2 Returning the Device and Parts.....	22
4.2 Selecting an Installation Site.....	10	12 Disposal	22
4.3 Unpacking.....	10	12.1 Decontaminating the Device.....	22
4.4 Acclimatizing.....	10	12.2 Disposing of the Device and Parts.....	22
5 Getting Started	10	13 Technical Data	23
5.1 Connecting the Power Supply.....	10	13.1 Dimensions and Weight.....	23
5.2 Connecting the Device to the Testing or Laboratory		13.2 Ambient Conditions at the Installation Site.....	23
Device (YDP30 Only).....	11	13.3 Power Supply.....	23
5.3 Connecting the Device to the Wi-Fi Router or		13.3.1 Device.....	23
Corporate Network (YDP30-NET Only).....	11	13.3.2 Power Supply Unit.....	23
5.4 Switching the Device On or Off.....	12	13.4 Operating Modes and Connection Types.....	24
5.5 Opening or Closing the Top Cover.....	12	13.5 Connections and Connection Adjustment Options	
5.6 Inserting Media into the Device.....	12	(YDP30 Only).....	24
5.6.1 Inserting Media into the Media Mount.....	12	13.6 Printing Process and Print Data.....	24
5.6.2 Inserting Release Paper of Self-adhesive		13.7 Ink Ribbon.....	25
Media (Peel-off Mode).....	13	13.8 Media.....	25
5.6.3 Setting the Print Roller to the		13.9 Acoustic Emission.....	25
Media Thickness.....	13		
5.6.4 Completing the Media Insertion.....	13		
5.7 Inserting the Ink Ribbon.....	14		
5.8 Printing Status Information (YDP30-NET Only).....	14		

14 Accessories, Consumables, and Spare Parts 25

14.1 Accessories 25

14.2 Consumables. 26

14.3 Spare Parts 26

15 Sartorius Service. 26

16 Conformity Documents 26

1 About this Document

1.1 Scope

These instructions are part of the device. These instructions apply to the device in the following versions:

Device	Type
YDP30, Standard printer YDP30-NET, Network printer	YDP30

1.2 Target Groups

These instructions are written for the following target groups. The target groups must have acquired the knowledge indicated below.

Target group	Knowledge and qualifications
User	The user is familiar with the operation of the device and the associated work processes. The user understands the hazards which may arise when working with the device and knows how to prevent them. They have been trained in the operation of the device.
Operator	The operator of the device is responsible for ensuring compliance with workplace health and safety regulations. The operator must ensure that all persons who work with the device have access to the relevant information and are trained in working with the device.

1.3 Symbols Used

1.3.1 Warnings in Operation Descriptions

WARNING

Denotes a danger with the risk that death or severe injury may result if it is **not** avoided.

NOTICE

Denotes a danger with the risk that property damage may result if it is **not** avoided.

1.3.2 Other Symbols

- ▶ Required action: Describes actions that must be carried out.
- ▷ Result: Describes the result of the actions carried out.

2 Safety Instructions

2.1 Intended Use

The device is intended to be used to print data from Sartorius testing or laboratory devices, for example for printouts of calibration or measurement results from integrity testing devices or laboratory balances. The device can be used in laboratories.

Printing is carried out using the "thermo transfer printing" or "thermo direct printing" printing processes. The device automatically selects the appropriate printing process based on the inserted media, for example standard paper or thermal paper.

Depending on the model, the device must be set up as either a standard printer or a network printer and connected to the testing or laboratory devices (suitability, see Chapter "13.4 Operating Modes and Connection Types", page 24).

The device is exclusively designed for use in accordance with these instructions. Any further use beyond this is considered **improper**.

If the device is **not** used properly: The protective systems of the device may be impaired. This can lead to unforeseeable personal injury or property damage.

Operating Conditions for the Device

Do **not** use the device in potentially explosive environments. Only use the device indoors.

The device may only be used with the equipment and under the operating conditions described in the Technical Data section of these instructions.

2.1.1 Modifications to the Device

If the device is modified, for example by attaching extra components: The safety of the device may be impaired or the device compliance may lose its validity.

If you have any queries regarding modifications to the device, contact Sartorius.

2.1.2 Device Repairs and Maintenance

Device repairs and maintenance may only be carried out by persons with appropriate specialized knowledge. If the device is **not** repaired or maintained by such a person: The safety of the device may be impaired or the quality marks may lose their validity.

Even outside the warranty period, we recommend that repairs be carried out by Sartorius Service or after consulting Sartorius Service.

Only the maintenance tasks described in these instructions should be carried out. For maintenance tasks that need to be carried out by Sartorius Service, contact Sartorius Service.

2.2 Personnel Qualification

If people who do **not** have sufficient knowledge on the safe handling of the device carry out work on the device: Those persons may injure themselves or others nearby.

- ▶ Ensure that all individuals working on the device possess the necessary knowledge and qualifications (description, see Chapter "1.2 Target Groups", page 5).
- ▶ If a particular qualification is required for the actions described: Allow these activities to be carried out by the required target group.
- ▶ If **no** particular qualification is required for the activities described: Have these tasks carried out by the "user" target group.

2.3 Significance of these Instructions

Failure to follow the instructions in this manual may have serious consequences, e.g. danger to individuals.

- ▶ Read the instructions carefully and in full.
- ▶ Ensure that the information contained in these instructions is available to all individuals working on the device.
- ▶ Retain the instructions.
- ▶ If these instructions are lost, request a replacement or download the latest version from the Sartorius website (www.sartorius.com).

2.4 Functionality of the Device

A damaged device or worn parts may lead to malfunctions or cause hazards which are difficult to recognize.

- ▶ Only operate the device when it is safe and in perfect working order.
- ▶ Comply with the maintenance intervals (for intervals and maintenance work, see Chapter "7.2 Maintenance Schedule", page 16).
- ▶ Allow any damage to be repaired immediately by Sartorius Service.

2.5 Electrical Equipment

2.5.1 Damage to the Device's Electrical Equipment

Damage to the device's electrical equipment, e.g. damaged insulation, can be life-threatening. There is a danger to life from contact with live parts.

- ▶ If the electrical equipment of the device is defective, cut off the power supply and contact Sartorius Service.
- ▶ Keep live parts away from moisture. Moisture can cause short-circuits.

2.5.2 Power Supply Unit and Power Supply Cable

Serious injury can result, e.g. from electric shocks, if an unsuitable/inadequately dimensioned power supply cable or unsuitable power supply unit is used.

- ▶ Only use the original power supply unit and original power supply cable.
- ▶ If the power supply unit or power supply cable needs to be replaced: Please contact Sartorius Service. Do **not** repair or modify the power supply unit or power supply cables.

2.6 Conduct in an Emergency

If an emergency occurs, e.g., due to the malfunctions of the device or dangerous situations: People might be injured. The device must be immediately taken out of operation:

- ▶ Disconnect the power supply cable from the power supply at the installation site.
- ▶ Secure the device to prevent it from restarting.

2.7 Accessories, Consumables, and Spare Parts

Unsuitable accessories, consumables, and spare parts can adversely affect functionality and safety and can have the following consequences:

- Risk of injury to persons
 - Damage to the device
 - Device malfunctions
 - Device failure
-
- ▶ Only use approved accessories, consumables, and spare parts supplied by Sartorius.
 - ▶ Only use accessories, consumables, and spare parts that are in proper working order.

3 Device Description

3.1 Device Overview

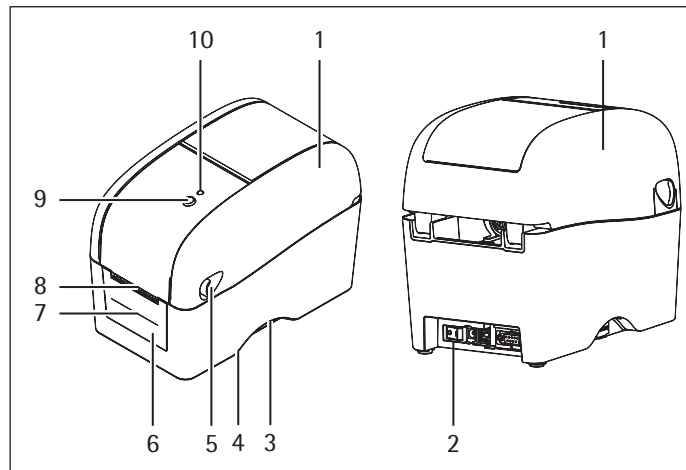
Fig. 1: YDP30 (front and rear view, example)

Pos.	Name	Description
1	Top cover	
2	On/Off switch	
3	Recessed grip	
4	Manufacturer's ID label	Attached to the underside of the device.
5	Lever	Is intended to open the top cover.
6	Front cover	
7	Release paper opening	For the output of release paper from self-adhesive media (peel-off mode).
8	Paper output tray	
9	Feed button	- Is intended to move inserted media forward, e.g. when printing labels. - Is intended to stop and resume printing. - Is intended to select certain functions, e.g. configuration of the RS232 interface, calibration of sensors.
10	LED display	

3.2 Connections

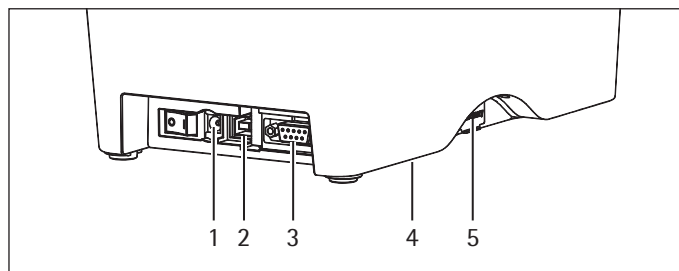
Fig. 2: Connections (rear view YDP30, example)

Pos.	Name	Description
1	Power connection	Is intended to connect the power supply cable.
2	USB	USB port
3	RS-232C	Serial interface, YDP30 only
	Ethernet	YDP30-NET only
4	Production interface	No function
5	Interface for SD memory card	No function

3.3 LED Display

The LED display denotes the device's operating status.

Light color and behavior	Description
Green	The device is ready for use.
Flashing green	The device is receiving data or a print job has been stopped.
Illuminated yellow/red	The device is connecting to the Wi-Fi network or corporate network (YDP30-NET only).
Illuminated yellow	The device is deleting internally stored data.
Illuminated red	The top cover is open.
Flashing red	There is a problem, e.g. a paper jam.

3.4 Components for Print Function and Lock for Top Cover

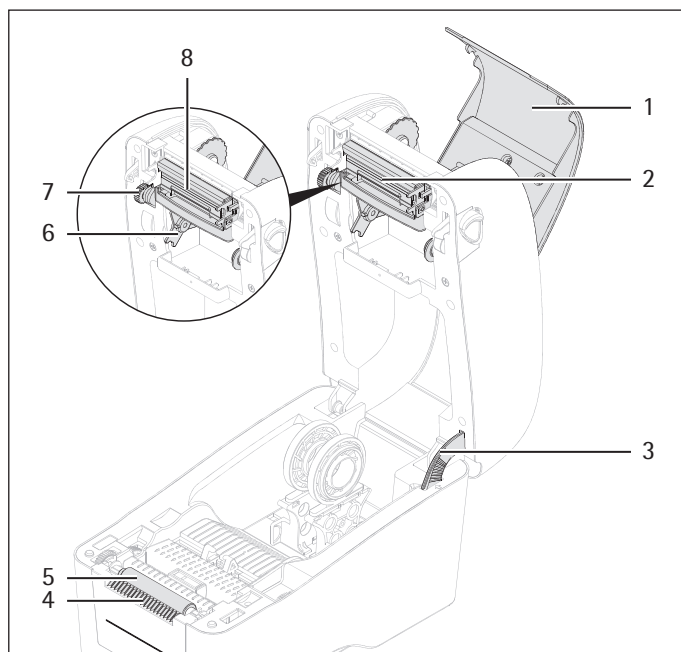
Fig.3: Components for print function and lock for top cover (top cover and ink ribbon cover open, example)

Pos.	Name	Description
1	Ink ribbon cover	
2	Print head	
3	Lock for top cover	Keeps the top cover in an open position. This means that the top cover cannot close accidentally.
4	Tear bar	
5	Print roller	
6	Intake for the ink ribbon supply spool	
7	Adjustment wheel to wind the ink ribbon	
8	Intake for the ink ribbon take-up spool	

3.5 Components for Media Feed and Sensors

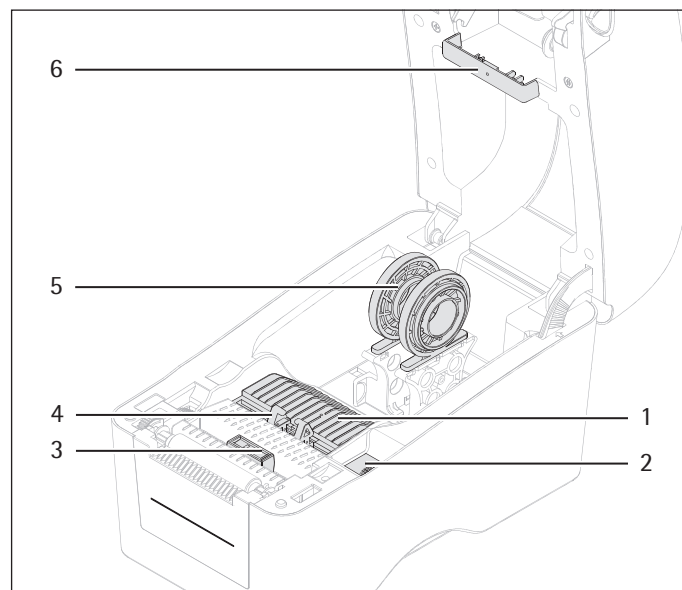
Fig. 4: Components for media feed (top cover open, example)

Pos.	Name	Description
1	Distance sensor (transmitter)	
2	Adjustment wheel for media feed	Sets the media feed to the width of the medium.
3	Sensor for black marking	
4	Media feed	Keeps the inserted media in position. Can be set to different widths.
5	Media mount	Holds the spool of the medium. Consists of 2 mounts, the distance between which can be changed. The mounts can be rotated 180°. This allows the media mount to be set to 2 media spool sizes.
6	Distance sensor (receiver)	

4 Installation

4.1 Scope of Delivery

Item	Quantity	Purpose
Device	1	
Power supply unit, type WDS060240	1	
Country-specific power supply cable YDP30 only	1	
Printer connection cable, type USB	1	
Printer connection cable, type RS232	1	
Winding roller	1	Used to wind the ink ribbon
Operating Instructions	1	

4.2 Selecting an Installation Site

Procedure

- Make sure that the following conditions are met at the installation site:

Condition	Features
Ambient conditions	Suitability tested (see Chapter "13.2 Ambient Conditions at the Installation Site", page 23)
Setup surface	– Stable, even surface with little vibration – Sufficient space for the device.

4.3 Unpacking

Procedure

- Unpack the device.
- Install the device at the intended installation site. The device must **not** be dropped or shaken.
- We recommend you keep the original packaging for a safe and proper return of the device, e.g. for repairs.

4.4 Acclimatizing

When a cold device is brought into a warm environment: The temperature difference can lead to condensation from humidity in the device (moisture formation). Moisture in the device can lead to malfunctions.

Procedure

- Allow the device to acclimatize for approx. 2 hours at the installation site. The device must be disconnected from the power supply during that time.

5 Getting Started

5.1 Connecting the Power Supply

Procedure

-  **WARNING** Severe injuries may be caused by using defective power supply cables! Check the power supply cable for damage, e.g., cracks in the insulation.
 - If required: Please contact Sartorius Service.
- Check whether the country-specific power plug matches the power connections at the installation site.
 - If required: Please contact Sartorius Service.
- Check whether the voltage specifications on the manufacturer's ID label match those of the power supply at the installation site.
 - If the input voltage is too high: Do **not** connect the device to the power supply.
 - Please contact Sartorius Service.
- Connect the power supply unit connection cable to the "Power Supply" connection.
- Connect the power supply cable to the power supply unit and the power supply at the installation site. The power outlet must be easy to reach.

5.2 Connecting the Device to the Testing or Laboratory Device (YDP30 Only)

Procedure

- Connect one of the printer connection cables to the appropriate port on the device and the testing or laboratory device (for information about ports on the testing or laboratory device, see the testing or laboratory device's operating instructions).

5.3 Connecting the Device to the Wi-Fi Router or Corporate Network (YDP30-NET Only)

Multiple testing and laboratory devices can access the device. Depending on the configuration, the testing or laboratory devices can access the device via an independent Wi-Fi router or via a corporate network.

Connection Type	Connection of Device and Testing or Laboratory Device	Connection Type Advantage	Connection Type Disadvantage
Access via Wi-Fi router (independent Wi-Fi network)	– The device is connected to the Wi-Fi router via a network connection cable. – The testing or laboratory device is connected to the Wi-Fi router via Wi-Fi components.	– Printing can take place via an independent Wi-Fi network. – The device can be used without considering the corporate network's security guidelines or network settings.	Only the Wi-Fi print function of the testing or laboratory device can be used.
Access via corporate network	– The device is connected to the corporate network via a network connection cable. – The testing or laboratory device is connected to the corporate network via Wi-Fi components.	– Fewer connection cables are required. – No Wi-Fi router is required. – All Wi-Fi functions of the testing or laboratory device can be used.	The device does not support fixed IP addresses. For this reason, the Wi-Fi functions of the testing or laboratory devices may be restricted, e.g. only standard authentication methods are available.

Requirements

- The device and the Wi-Fi router are turned off.
- The testing or laboratory device is enabled and prepared for Wi-Fi operation (see the testing or laboratory device's operating instructions).

Procedure

- Check the suitability of the Wi-Fi router or the corporate network (for suitability see Chapter "13.4 Operating Modes and Connection Types", page 24).
- Connect a network connection cable to the "Ethernet" port and to the Wi-Fi router or corporate network.

5.4 Switching the Device On or Off

Procedure


- ▶ To switch the device on: Press the On/Off switch so that it is in the "I" position.
- ▷ When the device is connected to a network: The device establishes a connection to the Wi-Fi network or corporate network and retrieves an IP address for the device. This procedure usually takes 20 seconds.
- ▷ The LED display illuminates in the color "yellow".
- ▷ If media and an ink ribbon are inserted into the device: The device will print an initial printout containing the Sartorius logo and the current software version.
- ▷ The LED display illuminates in the color "green". The device is ready for use.
- ▶ To switch the device off: Press the On/Off switch so that it is in the "O" position.

5.5 Opening or Closing the Top Cover

Requirements

The device is switched off.

Procedure

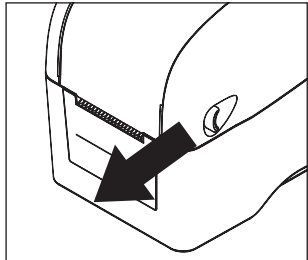- ▶ To open the top cover:
 - ▶ Pull the levers located on both sides of the top cover forwards.

- ▶ Open the top cover upward and push it backward until it stops.
- ▷ The top cover lock audibly engages.
- ▶ To close the top cover:
 - ▶ Push sideways against top cover lock. This will unlock the top cover.
 - ▶ Push the top cover back down. Both levers should snap into place.

5.6 Inserting Media into the Device

5.6.1 Inserting Media into the Media Mount


The media spool, e.g. the standard paper roll, must be inserted between the 2 media mounts. The print side of the media must be facing up.

Requirements

- The print roller is clean (see Chapter "7.2 Maintenance Schedule", page 16).
- The device is switched off.


Procedure

- ▶ Check if the media are suitable for the device (suitability, see Chapter "13.7 Ink Ribbon", page 25).
- ▶ Open the top cover.

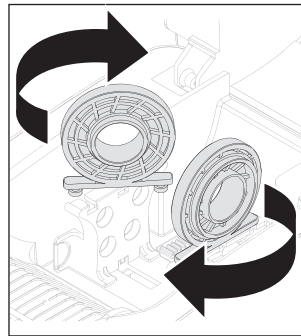- ▶ Check if the media mount is set to the right media spool size.
 - ▶ If required: Push the 2 mounts apart and turn them 180°.

- ▶ Push the 2 mounts apart and fix the media spool between the mounts.

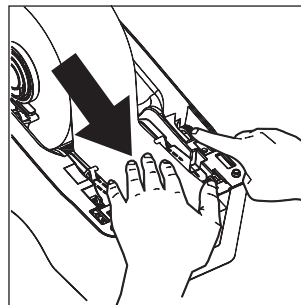- ▶ Push the media over the media feed, the sensor for black marking and the print roller.

- ▶ Turn the setting wheel for the media feed so that the media feed is set to the width of the inserted media.
- ▶ Close the top cover.

5.6.2 Inserting Release Paper of Self-adhesive Media (Peel-off Mode)

If the release paper of self-adhesive media has been inserted: The release paper will automatically be removed during printing.

Requirements

The device is switched off.

Procedure

- Open the top cover.

- Open the front cover.
- Separate the release paper from the media and feed it through the release paper opening.

- Close the front cover. The front cover must snap into place.
- Close the top cover.

5.6.3 Setting the Print Roller to the Media Thickness

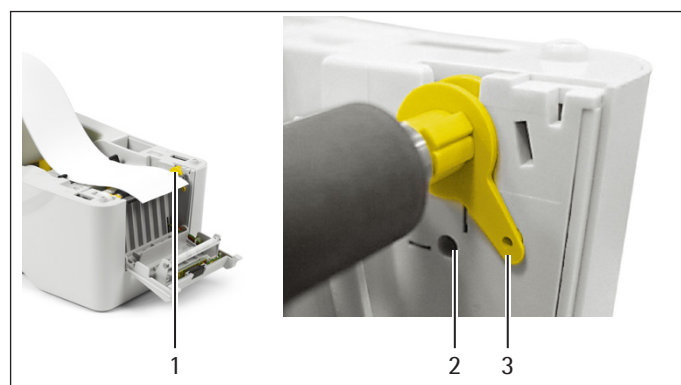
Fig. 5: Lever (set in the position for normal paper, example)

Pos.	Name	Description
1	Lever	For setting the print roller.
2	Position for normal paper	For setting the print roller to normal paper thicknesses (for media thickness, see Chapter "13.8 Media", page 25).
3	Position for thick paper	For setting the print roller to the maximum paper thickness (for maximum paper thickness, see Chapter "13.8 Media", page 25).

Requirements

The device is switched off.

Procedure

- Open the top cover.
- Fold the front cover downward.
- Push the lever to the desired position, e.g. into the position for thick paper.
- Fold the front cover upward and let it snap into place.
- Close the top cover.

5.6.4 Completing the Media Insertion

Procedure

- Define the media sensor and calibrate the sensors (see Chapter "7.6 Calibrating Sensors and Changing Settings", page 18).


5.7 Inserting the Ink Ribbon

Requirements


- The print head is clean (see Chapter "7.2 Maintenance Schedule", page 16).
- The device is switched off.

Procedure


- ▶ Open the top cover.
- ▶ Open the ink ribbon cover.

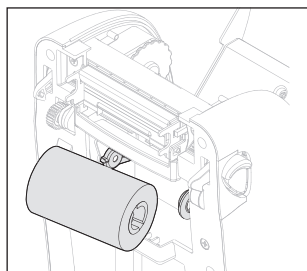- ▶ Insert the ink ribbon supply spool into the intake for the ink ribbon supply spool starting with the right side, then the left side, and then align it.

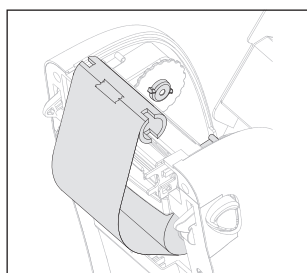- ▶ Attach the ink ribbon to the supplied winding roller, e.g. with tape.
- ▶ Pull the ink ribbon through the winding roller over the print head and place first the right side and then the left side of the winding roller on the winding spool for the ink ribbon and align it with the recesses.

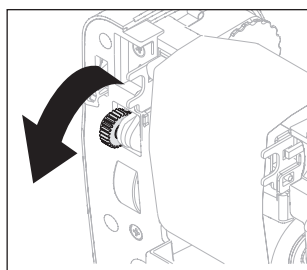- ▶ Turn the setting wheel downward to wind the ink ribbon until the plastic leader of the ink ribbon is tightly wound and the black part of the ink ribbon is on the print head.

- ▶ Close the ink ribbon cover.
- ▶ Close the top cover.

5.8 Printing Status Information (YDP30-NET Only)

The device's status information is needed for configuring the testing or laboratory device. The status information should print within 10 seconds of turning on the device. The printout should contain the following information:

- Line "IP ADDRESS": IP address retrieved from the Wi-Fi network or corporate network by the device.
- Line "HOST NAME": Host name retrieved from the Wi-Fi network or corporate network by the device.

If the status information **cannot** be printed: There is a fault in the device or Wi-Fi router, or the connection to the corporate network is interrupted.

Requirements

If a Wi-Fi router is used: The Wi-Fi router is turned on and ready for use.

Procedure

- ▶ If the LED display on the device illuminates in green: Switch the device off and on again.
- ▶ If the status printout has printed and while the LED display is illuminated in "yellow": Press the feed button.
- ▷ The device will print the device status information. The IP address of the network printer is specified in the "IP ADDRESS" line. The host name is specified in the "HOST NAME" line.
- ▶ If the IP address entry is "0.0.0.0":
 - ▶ Press the feed button again.
 - ▷ The status information will be printed again.
 - ▷ If required: Repeat this procedure up to 5 times.
- ▷ If the status information **cannot** be printed or the IP address or host name information are still missing:
 - ▶ Check and adjust the settings for the Wi-Fi server or corporate network.
 - ▶ If required: Please contact Sartorius Service.
- ▶ Keep the printout of the status information for the configuration of the testing or laboratory device.

6 System Settings

6.1 Configure "RS232" Connection (YDP30 Only)

The "RS232" connection can be configured in 2 configuration options (see Chapter "13.5 Connections and Connection Adjustment Options (YDP30 Only)", page 24). The feed button is used to carry out the configuration.

Requirements

The device is ready for use.

Procedure

- ▶ Switch off the device.
- ▶ Switch on the device again.
- ▶ Once the initial printout has printed: Press the feed button within the next 3 seconds. This starts the configuration mode.
- ▷ A configuration routine is started, the configuration options are printed. The following data is printed:

Configuration options	Description
SETUP SERIAL PORT:	
PRESS KEY: XXXXBdYZZ	Shows the adjustable configuration, e.g. 9600BD801. The configuration can be selected by pressing the feed button.
WAIT 5s: no change	Wait 5 seconds for the next option to be displayed.

- ▶ If you want to select the printed configuration option: Press the feed button within the next 5 seconds.
- ▷ If the configuration option is **not** selected: The next configuration option will be printed in the configuration routine.
- ▷ If the configuration option is **not** selected again: The configuration routine will be canceled: The notification [no changes made/EXIT SETUP/print settings] will be printed.
- ▷ The device's configuration data will be printed (see Chapter "7.6.1 Check the Configuration Data of the Device", page 19).

7 Cleaning and Maintenance

7.1 Cleaning Device Surfaces

Procedure

- ▶ Switch off the device.
- ▶ Disconnect the power supply cable from the power supply at the installation site.
- ▶ **NOTICE** Corrosion or damage to the device due to unsuitable cleaning agents!
 - ▶ Do **not** use corrosive, chloride-containing or aggressive cleaning agents.
 - ▶ Do **not** use cleaning agents that contain abrasive ingredients, e.g. scouring agents, steel wool.
 - ▶ Do **not** use solvent-based cleaning agents.
 - ▶ Make sure that cleaning materials, such as cloths, are only slightly damp. Do **not** use wet cleaning cloths.
- ▶ Wipe the housing with a slightly damp cloth. Use a mild soapy solution or a suitable cleaning agent for more severe contaminations.

7.2 Maintenance Schedule

Interval	Component	Action	Chapter, page
Before media are inserted	Print head	Remove the ink ribbon.	7.3, 16
		Clean the print head.	7.5.1, 17
	Print roller	Remove the media.	7.4, 17
		Clean the print roller.	7.5.2, 17
After media are inserted	Sensors	Calibrate the sensors.	7.6, 18
Monthly	Sensors	Remove the media.	7.4, 17
		Clean the sensors.	7.5.3, 17
As required	Tear bar	Remove the media.	7.4, 17
		Clean the tear bar.	7.5.2, 17
	Interior	Remove the media.	7.4, 17
		Clean the interior components.	7.6, 18

7.3 Removing the Ink Ribbon

Procedure

- ▶ Switch off the device.
- ▶ Open the top cover.
- ▶ Open the ink ribbon cover.
- ▶ Remove the ink ribbon's paper winding spool from the winding spool for the ink ribbon first on the left side, then the right.
- ▶ Remove the ink ribbon's supply spool from the intake for the ink ribbon's supply spool first on the left side, then the right.

7.4 Removing Media from Media Mounts

Procedure

- Switch off the device.
- Open the top cover.
- If release paper of self-adhesive media has been inserted:
 - Open the front cover.
 - Pull the release paper out of the release paper opening.
 - Close the front cover.
- Push the 2 mounts of the media feed apart and remove the media spool.

7.5 Cleaning the Sensors and Interior Components

7.5.1 Cleaning the Print Head

The print head must **not** be touched. The print head must be cleaned with 100% ethanol and a cotton swab or with a print head cleaning pen.


The print head must **not** come into contact with medical alcohol. Medical alcohol can damage the print head.

Requirements

An ink ribbon has **not** been inserted into the device.

Procedure

- Switch off the device.
- Open the top cover.
- Open the ink ribbon cover.

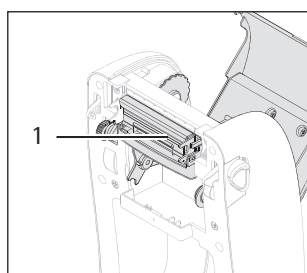- Use a cotton swab moistened with 100% ethanol or a print head cleaning pen to clean the print head (1).

- If the print head accidentally comes into contact with skin: Clean the print head again.

7.5.2 Cleaning the Print Roller or Tear Bar


The print roller and tear bar must be cleaned with 100% ethanol and a cotton swab or lint-free cloth.

Requirements

There are **no** media in the device.

Procedure

- Switch off the device.
- Open the top cover.
- Moisten a cotton swab or lint-free cloth with 100% ethanol.

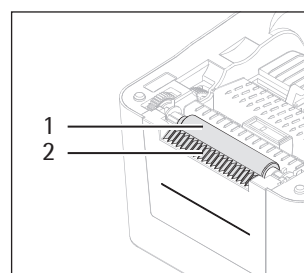- Wipe the print roller (1). Turn the print roller so that the entire surface of the print roller can be cleaned.
- Wipe the tear bar (2).

7.5.3 Cleaning Sensors or the Interior


The sensors and the interior of the device must be cleaned with an air brush or a vacuum brush.

Requirements

There are **no** media in the device.

Procedure

- Switch off the device.
- Open the top cover.

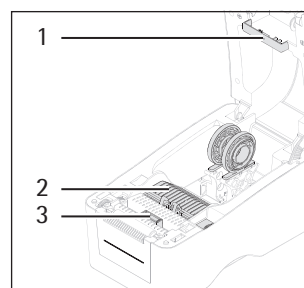- Clean the distance sensor (receiver) (1), distance sensor (transmitter) (2), and sensor for black marking (3) with an air brush or a vacuum brush.

- Clean the interior components with an air brush or a vacuum brush.

7.6 Calibrating Sensors and Changing Settings

The calibration of sensors and setting of the device functions is done using the feed button while the device is being switched on. The LED display will flash and change colors throughout the process.

Color sequence (change every 5 seconds)

Red	Yellow	Green	Green/yellow	Red/yellow
-----	--------	-------	--------------	------------

The feed button must be held down while the colors change until the LED display shows the correct color for the desired function, e.g. "red" to calibrate the sensors.

Available settings	LED display color	Description
Define the sensor for black marking as the media sensor and calibrate	Green/yellow	To calibrate labels with black markings
Define the distance sensor as the media sensor and calibrate	Red/yellow	To calibrate labels with gaps
Calibrate sensor for black marking and distance sensor	Red	For the initial calibration when using individual labels. The calibration is saved to protected memory.
Calibrate sensor for black marking and distance sensor, perform self-test, and print in core dump mode	Yellow	The following functions are run one after the other: – The sensitivity of the sensors is set. During this process, the device measures the length of the labels. – The device performs a self-test and prints the configuration data (see Chapter "7.6.1 Check the Configuration Data of the Device", page 19). – The device prints a core dump (see Chapter "7.6.2 Checking Core Dump", page 19).
Restore device to factory settings	Green	If the device has been calibrated to print individual labels: To change calibration to print on continuous paper The following functions are conducted: – The internal memory (DRAM) is deleted. – The device is reset to factory settings.

Requirements

- The ink ribbon and media have been inserted in the device.
- If a self-test needs to be performed: The inserted media have the required minimum width (see Chapter "12.1 Decontaminating the Device", page 22).

Procedure

- Switch off the device.
- Press and hold the feed button while switching on the device.
- Hold the feed button down until the LED display shows the correct color for the desired function, e.g. "red" to calibrate the sensors.
- When the LED display shows the correct color: Release the feed button.

- ▷ The device will perform the functions, e.g. calibration of sensors.
- ▷ Depending on the device function, media will come out of the paper output tray with or without print function.

7.6.1 Check the Configuration Data of the Device

If a self-test has been performed or settings have been changed: The device's configuration data will be printed. The printout contains the following information:

Configuration data	Description
YDP30 YDP30-NET Version: X.XX. XX	Name and firmware version on the main board
SERIAL NO.:	No entry
MILAGE (m): XXX	Already printed (in meters)
CHECKSUM: XXXXXXXX	Checksum of the main board firmware
SERIAL PORT: XXXX, X , X , X	Configuration of the "RS232" connection
CODE PAGE: XXX	Code page
COUNTRY CODE: XXX	Country code
SPEED: X INCH	Print speed in inch/second
DENSITY: X. X	Print density
SIZE: X.XX , X.XX	Label size: Width, height
GAP: X.XX , X.XX	Gap size (vertical gap, offset)
TRANSPARENCE: X	Sensor sensitivity
[FILE LIST:]	File management information, e.g. number of files in the internal memory
[Sample]	Print head test sample

Procedure

- ▶ Check the information on the printout.
 - ▶ If required: Change the settings.

7.6.2 Checking Core Dump

In the core dump, characters are listed in two columns. The received characters are listed in the left column. The hexadecimal values of the characters are listed in the right column.

Procedure

- ▶ Check the information on the printout.
 - ▶ If required: Change the settings.

7.7 Re-starting the System

Procedure

- ▶ Commission the device (see Chapter "5 Getting Started", page 10).

8 Malfunctions

8.1 Device Malfunctions

Malfunction	Cause	Solution	Chapter, page
The LED display is inactive.	The device is switched off.	Switch on the device.	5.4, 12
	The device is not connected to the power supply at the installation site.	Connect the device to the power supply at the installation site.	5.1, 10

8.2 Malfunctions when Printing

Malfunction	Cause	Solution	Chapter, page
The device does not print.	A print job has been stopped. The LED display flashes green.	Press the feed button. This will resume the print job.	
	No media have been inserted. The LED display flashes red.	Insert the media.	5.6.1, 12
	The inserted media have been used up.	Remove the media spool.	7.4, 17
		Insert the media.	5.6, 12
		Press the feed button. This will resume the print job.	
	The ink ribbon has been inserted incorrectly.	Remove the ink ribbon.	7.3, 16
		Insert the ink ribbon.	5.7, 14
The device does not print (only when used as standard printer).	The device is not correctly connected to the testing or laboratory device.	Check the connection.	5.2, 11
	The RS232 connection cable is not correctly assigned.	Use the supplied printer connection cable, type RS232.	
		If the printer connection cable is missing: Please contact Sartorius Service.	
	The settings for the serial interface between the device and the testing or laboratory device differ.	Change the configuration for the "RS232" connection.	6.1, 15
The device does not print (only when used as network printer).	The device received multiple print jobs simultaneously and has only printed the first print job.	Restart the missing print job.	
	The device is not correctly connected to the testing or laboratory device.	Print the status information.	5.8, 14
		Check if the device is correctly configured on the testing or laboratory device (see the testing or laboratory device's operating instructions).	
The media shift during printing.	The printer settings are incorrect.	Change the media sensor settings.	7.6, 18
		Reset the device to factory settings.	7.6, 18
		Calibrate the sensors.	7.6, 18

Malfunction	Cause	Solution	Chapter, page
There is a paper jam. The LED display flashes red.	The sensors are incorrectly set.	Calibrate the sensors.	7.6, 18
	Media are jammed between components.	Open the top cover.	
		Remove the jammed media.	

8.3 Poor Print Quality

Malfunction	Cause	Solution	Chapter, page
The print quality is poor.	The top cover is not completely closed.	Close the top cover. Both levers should snap into place.	
	The media have not been inserted correctly.	Remove the media.	7.4, 17
		Insert the media.	5.6, 12
	The print head is dirty.	Clean the print head.	7.5.1, 17

8.4 Re-starting the System

Procedure

- Commission the device (see Chapter "5 Getting Started", page 10).

9 Decommissioning

9.1 Decommissioning the Device

Procedure

- Switch off the device.
- Disconnect the device from the power supply and all devices. In order to do this, remove all connection cables.
- Remove media from media mounts.
- Remove the ink ribbon.
- Clean the device.

10 Transport

10.1 Transporting the Device

Procedure

- ▶ Switch off the device.
- ▶ Disconnect the device from the power supply and all devices. In order to do this, remove all connection cables.
- ▶ Grab the device by the recessed grips and carry to the desired installation site.

11 Storage and Shipping

11.1 Storing

Requirements

The device has been decommissioned.

Procedure

- ▶ Store the device according to the ambient conditions (see Chapter 7.6, page 18).

11.2 Returning the Device and Parts

Defective devices or parts can be sent back to Sartorius. Returned devices must be clean, decontaminated, and properly packed.

Transport damage as well as measures for subsequent cleaning and disinfection of the device or parts by Sartorius shall be charged to sender.

Devices contaminated with hazardous materials, e.g. harmful biological or chemical substances, will **not** be accepted for repair or disposal. The devices must be decontaminated before shipping (for decontamination, see Chapter "12.1 Decontaminating the Device", page 22).

Procedure

- ▶ Decommission the device.
- ▶ Contact Sartorius Service for instructions on how to return devices or parts (please refer to www.sartorius.com).
- ▶ Pack the device and its parts properly for return.

12 Disposal

12.1 Decontaminating the Device

The device does **not** contain any hazardous materials that necessitate special disposal measures. If the device has come into contact with hazardous substances: Steps must be taken to ensure proper decontamination and declaration.

Procedure

- ▶ If the device has come into contact with hazardous substances: Decontaminate the device. The operator of the device is responsible for adhering to local government regulations on the proper decontamination and declaration for transport and disposal.

12.2 Disposing of the Device and Parts

The device and the device accessories must be disposed of properly by disposal facilities.

A lithium battery is installed inside the device. Lithium batteries must be disposed of properly by disposal facilities.

The packaging is made of environmentally friendly materials that are intended to be used as secondary raw materials.

Requirements

The device has been decontaminated.

Procedure

- ▶ Dispose of the device. Follow the disposal instructions on our website (www.sartorius.com). Inform the disposal facility that there is a lithium battery installed inside the device.
- ▶ Dispose of the packaging in accordance with local government regulations.

13 Technical Data

13.1 Dimensions and Weight

	Unit	Value
Dimensions (L x W x H)	mm	241.3 x 139.9 x 177.4
Weight	kg	2.2

13.2 Ambient Conditions at the Installation Site

	Unit	Value
Conventional laboratory rooms, max. height above sea level	m	2000
Temperature		
Operation	°C	+5 – +40
Storage	°C	-40 – +60
Relative humidity, non-condensing		
Operation	%	25 – 85
Storage	%	10 – 90

13.3 Power Supply

13.3.1 Device

	Unit	Value
Voltage	V _{DC}	24
Current consumption	A	2.5
Power supply only by Sartorius power supply unit WDS060240		

13.3.2 Power Supply Unit

	Unit	Value
Type		
Sartorius power supply unit only, model WDS060240		
Primary		
Voltage	V _{AC}	100 – 240
Frequency	Hz	50 – 60
Current consumption, maximum	A	1.6
Secondary		
Voltage	V _{DC}	24
Current, maximum	A	2.5

13.4 Operating Modes and Connection Types

Use as standard printer

YDP30

Connection to the testing or laboratory device via printer connection cable

Use as network printer

YDP30-NET

Connection to the Wi-Fi router or corporate network via network connection cable

Wi-Fi router requirements: Hotspot function and DHCP server must be available

Corporate network requirements: DHCP server must be present

13.5 Connections and Connection Adjustment Options (YDP30 Only)

	Unit	Value
Speed of connections		
"USB" connection, type 2.0, full speed		
"RS232" connection, maximum	bps	115,200
Standard setting for "RS232" connection		
Setting: 9600BD801/9600.0.8.1		
Baud rate 9600/parity: Odd/8 data bits/1 stop bit		
Alternate setting for "RS232" connection		
Setting: 1200BD701/1200.0.7.1		
Baud rate 1200/parity: Odd/7 data bits/1 stop bit		

13.6 Printing Process and Print Data

	Unit	Value
Possible print process, according to the inserted media		
Thermo transfer printing		
Thermo direct printing		
Resolution	dpi	203
Print width, maximum	mm	54

13.7 Ink Ribbon

	Unit	Value
Ink ribbon output, maximum	m	90
Outer diameter of the ink ribbon, maximum	mm	37
Outer diameter of the ink ribbon spool, maximum	mm	12.7
Ink ribbon spool width	mm	56 – 58

13.8 Media

	Unit	Value
Suitability		
Standard paper		
Thermal paper, standard or self-adhesive		
Labels, size 58 mm x 30/76/100 mm		
Media with gaps or black markings		
Type of winding: Outer winding		
Dimensions		
Outer diameter of the media, maximum	mm	127
Media width	mm	15 – 60
Media thickness	mm	0.06 – 0.19
Outer diameter of the media spool	mm	25 – 38
Minimum width of the media for the self-test	mm	5.08

13.9 Acoustic Emission

	Unit	Value
Acoustic pressure level, maximum	dB (A)	<70

14 Accessories, Consumables, and Spare Parts

14.1 Accessories

This table contains an excerpt of the accessories that can be ordered. For information on other products, contact Sartorius.

Item	Quantity	Order number
Wireless nano router, for use as an independent Wi-Fi network	1	YWLAN02MS

14.2 Consumables

This table contains an excerpt of the consumables that can be ordered. For information on other products, contact Sartorius.

Item	Quantity	Order number
Set: Standard paper and ink ribbon	1	69Y03285
Set: Self-adhesive paper and ink ribbon	1	69Y03286
Standard thermal paper	1	69Y03287
Self-adhesive thermal paper	1	69Y03288
Self-adhesive labels		
58 mm x 100 mm, 350 pieces	1	69Y03094
58 mm x 76 mm, 500 pieces	1	69Y03093
58 mm x 30 mm, 1000 pieces	1	69Y03092

14.3 Spare Parts

This table contains an excerpt of the accessories that can be ordered. For information on other products, contact Sartorius.

Item	Quantity	Order number
Power supply unit WDS060240	1	69Y03294
Country-specific power supply cable		
EU	1	6900900
US/CDN	1	6900901
UK	1	6971945
AUS	1	6900905
ZA	1	6900902
Israel	1	6971975
Denmark	1	6971980
Italy	1	6971776
Printer connection cable, type mini USB to USB B	1	69Y03293
Printer connection cable, type RS232	1	69Y03295

15 Sartorius Service

Sartorius Service is always available to answer any queries regarding the device. Please visit the Sartorius website (www.sartorius.com) for information about service addresses, services provided, or to contact a local representative.

When contacting Sartorius Service with questions about the system or in the event of malfunctions, be sure to have the device information, e.g. serial number, close at hand. This information can be found on the manufacturer's type plate.

16 Conformity Documents

The attached documents confirm compliance of the device with the directives or standards cited.

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SARTORIUS

Manufacturer Sartorius Lab Instruments GmbH & Co. KG
37070 Goettingen, Germany

declares under sole responsibility that the equipment

Device type Printer

Model YDP30

in the form as delivered fulfils all the relevant provisions of the following British Regulations and meets the applicable requirements of the British Designated Standards including any amendments valid at the time this declaration was signed listed below:

The Electromagnetic Compatibility Regulations 2016
UK Statutory Instruments 2016 No. 1091
BS EN 61326-1:2013

The Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012
UK Statutory Instruments 2012 No. 3032
BS EN IEC 63000:2018

The Supply of Machinery (Safety) Regulations 2008
UK Statutory Instruments 2008 No. 1597
BS EN ISO 12100:2010, BS EN IEC 62368-1:2020*

The person authorised to compile the technical file:

Sartorius UK Ltd.
Longmead Business Centre, Blenheim Road
KT19 9 QQ Epsom, Surrey, UK

Sartorius Lab Instruments GmbH & Co. KG
Goettingen, 2021-08-16

Dr. Reinhard Baumfalk
Head of Product Development (LPS)

Halil Yildirim
Product Compliance Officer (SLI)

*: applied standard, which however is not harmonized for machines

FCC Supplier's Declaration of Conformity

SARTORIUS

Device type Printer

Model YDP30

Party issuing Supplier's Declaration of Conformity / Responsible Party – U.S. Contact Information

Sartorius Corporation
5 Orville Dr Suite 200
11716 Bohemia, NY, USA
Telephone: +1.631.254.4249

FCC Compliance Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Information to the user

Note: This equipment has been tested and found to comply with the limits for a **class B** digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Connections between the device and peripherals must be made using shielded cables in order to maintain compliance with FCC radio frequency emission limits.

Any modifications made to this device that are not approved by Sartorius may void the authority granted to the user by the FCC to operate this equipment.

Issued: 2021-08-25

Sartorius Lab Instruments GmbH & Co. KG
Otto-Brenner-Strasse 20
37079 Goettingen, Germany

Phone: +49 551 308 0
www.sartorius.com

The information and figures contained in these instructions correspond to the version date specified below.

Sartorius reserves the right to make changes to the technology, features, specifications and design of the equipment without notice.

Masculine or feminine forms are used to facilitate legibility in these instructions and always simultaneously denote all genders.

Copyright notice:

These instructions, including all components, are protected by copyright.

Any use beyond the limits of the copyright law is not permitted without our approval.

This applies in particular to reprinting, translation and editing irrespective of the type of media used.

Last updated:

10 | 2021